Quick Protocol



ScreenFect®siRNA Transfection Reagent

Package Contents

Cat. No.	ScreenFect®siRNA	Dilution Buffer
S-4001-2	0.2 ml	10 ml
S-4001	1.0 ml	50 ml
S-4001-3	5 x 1.0 ml	5 x 50 ml

Storage Conditions


Store ScreenFect® Reagents at 4°C. Do not freeze. For optimal long term activity, do not allow ScreenFect® Reagents to warm to room temperature each time it is used.

After several months of storage without using the reagent a slight precipitation might occur. If vortexed thoroughly, this has no influence on the performance of ScreenFect® Reagents.

General Considerations

For optimal results, amounts of ScreenFect®siRNA and nucleic acid (NA) need to be optimized for each cell type and each NA used. An optimization protocol is provided in our product manual which can be downloaded from the Incella homepage. We strongly recommend the One-Step transfection method for all of our products. For transfection of adherent cells, remove the used medium and mix fresh medium with the transfection complexes. Then add the mix to the cells.

ScreenFect® Protocol: One-Step Transfection

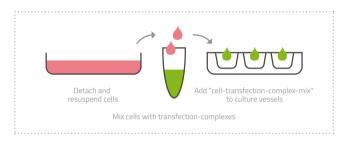

ScreenFect® Products

ScreenFect®A	Multipurpose reagent (most suitable for pDNA transfection, suitable for RNA applications) with very low cytotoxicity.	
ScreenFect®A-plus	Multipurpose reagent with optimized formulation requiring less reagent per transfection.	
ScreenFect®siRNA	Specialized reagent for best performance in siRNA delivery.	
ScreenFect®65	Reagent kit for protein production in HEK suspension cells.	

For additional information regarding ScreenFect®siRNA and other ScreenFect® Products, visit our homepage (www.screenfect.com) and view our product pages and instruction manuals.

Quick Protocol



Protocol for siRNA Transfection

	Component	Procedure for one well (96-well-plate)	96-well	24-well	6-well		
	Reagent Dilution	Dilute 0.25 µl of ScreenFect®siRNA in Dilution Buffer to a final volume of 10 µl and mix thoroughly.	0.25 μl reagent 10 μl dilution	1 μl reagent 40 μl dilution	4 μl reagent 120 μl dilution		
		Important: Vortex the reagent once per day of use. Add ScreenFect®siRNA reagent directly into supplied buffer with rapid pipette mixing or vortexing.					
2	siRNA Dilution	Dilute a total of 1 pmol siRNA in Dilution Buffer to a final volume of 10 $\mu l. $	1 pmol siRNA 10 µl dilution	5 pmol siRNA 40 µl dilution	25 pmol siRNA 120 µl dilution		
		Tip: Some siRNAs may require lower or higher amounts depending on the degree of expression of targeted gene.					
3	Complex formation	Combine the diluted siRNA with the ScreenFect®siRNA dilution and mix immediately using 10 rapid pipette strokes. Leave for 20 min at room temperature for complex formation.	20 µl complexes	80 µl complexes	240 µl complexes		
		Important: Do not vortex!					
4	Cell preparation & transfection	Add 80 μl freshly detached and resuspended cells to the complexes and mix with pipette.	Add 80 µl cell suspension	Add 420 µl cell suspension	Add 1250 µl cell suspension		
		Tip: The time-saving reverse cell transfection method may not be suited for all cell types. To transfect adherent cells, first remove and discard medium from cells, then add 80 µl fresh culture medium to transfection complexes, mix with pipette and immediately apply to cells.					
5	Cell plating	Transfer the cells and complexes to one well of a 96-well plate.	Transfer cells with complexes to plate	Transfer cells with complexes to plate	Transfer cells with complexes to plate		

Note: This protocol is a guideline. Values are suitable for easy for transfect cell lines. This protocol does not replace optimization experiments. View our manual for instructions. Serum does not affect the performance of ScreenFect®siRNA but we recommend avoiding antibiotics.

Cells must be mycoplasma free, in exponential growth phase and have even plating density across the entire surface area.